
ELECTRODO INFANTE MINDRAY MR61

MINDRAY MR61 INFANT ELECTRODE

SKU: 252303552028

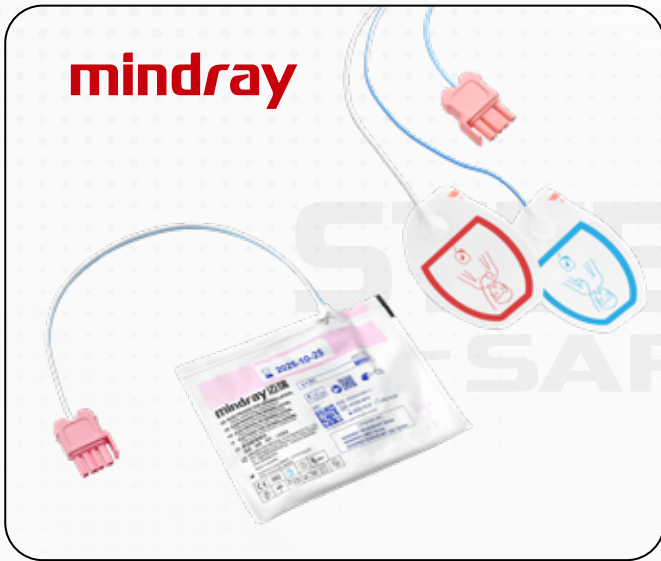DESCRIPCIÓN

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DESCRIPTION

The MR61 multifunction defibrillation electrodes are pre-gelled, self-adhesive electrodes designed for ECG monitoring, non-invasive pacing, manual defibrillation, and synchronized cardioversion in patients weighing less than 25 kg. This product must be used only by qualified clinical professionals and with specific Mindray defibrillators/monitors.

CHARACTERISTICS

The electrodes allow safe defibrillation with a maximum energy of up to 100 J. They must be kept in their original packaging until use, checking package integrity, gel condition, and expiration date, with a service life of 3 years. They are for single-patient, single-use only and must be replaced after 1 hour of continuous stimulation, 25 shocks, or every 8 hours. They must not be cleaned, sterilized, immersed in liquids, or applied to damaged skin, ensuring safe and reliable use. They incorporate a cable with a length of 120 cm.

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WARRANTY

If there is any manufacturing defect or inconformity you can contact the nearest distributor or write us directly to: asistencia@steelprosafety.com. The distributor is not liable to any injury, grievance or personal or patrimonial detriment resulting from improper use of the product. Before using this device make sure it is suitable, adequate or safe for the use intended.

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USES

They are designed for exclusive use with compatible Mindray equipment. Ideal for hospitals, schools, ambulances, and any environment that provides care for pediatric patients at cardiac risk.

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RESISTANCES

Before use, check the integrity of the packaging, the condition of the gel, and the expiration date. Do not crush, bend, or fold the electrodes, and do not use them if they show any damage, contamination, or dried gel. Do not open the package until application, do not apply additional gel, and do not clean, sterilize, or immerse them in liquids.